

Laboratory Approved Specimen Collection Devices & Media

For Corewell Health Lakeland Hospitals St. Joseph Laboratory,
Corewell Health Lakeland Hospitals Niles Laboratory and Corewell Health Lakeland Hospitals Watervliet Laboratory

Photos are for reference only. Always check container label, do not go by color alone!
General Test Listing. Not all tests are listed. Please review test catalog for more information.**



**Alinity m
multi-Collect Kit**
ITM-1190747

LAB1230885 Chlamydia+Gonococcus PCR✓
LAB1230890 STI Panel PCR** (CT, NG, TV, MG)✓
LAB1230886 Chlamydia PCR (CT)✓
LAB1230887 Gonococcus PCR (NG)✓

LAB1230889 Mycoplasma Genitalium PCR (MG)✓
LAB1230888 Trichomonas PCR (TV)✓

- ✓ **Sources for CT, NG, TV:** vaginal swabs (clinician-collected and self-collected in a clinical setting), endocervical swabs, gynecological specimens in ThinPrep® PreservCyt Solution, female urine, male urine and urethral swabs, oropharyngeal swabs, and rectal swabs on patients 14 years or older. Patients 13 years and younger must be collected using the APTIMA collection device. This testing will be sent out to a reference laboratory.
- ✓ **Sources for MG:** vaginal swabs (clinician-collected and self-collected in a clinical setting), endocervical swabs, male urine and urethral swabs on patients 14 years or older.



**ThinPrep® Pap
Container**
Not available via
Workday.
Order via
Client Services.

LAB1230097 Pap Test** (Indicate HPV Reflex**)
LAB1231770 HPV High Risk Screen with
HPV Genotype Reflex
LAB1230885 Chlamydia+Gonococcus PCR*

LAB1230886 Chlamydia PCR*
LAB1230887 Gonococcus PCR*
LAB1230888 Trichomonas PCR*

- ✓ Aseptic Control PAP Test Lubricating Jelly (ITM-1140484) and Surgilube® (ITM-1092842) are the only lubricating jellies to be used with ThinPrep®
- ✓ *Testing cannot be added onto a previous completed Pap Test Order. Note STI Panel and Mgen are not available for ThinPrep® Collection



**BD ESwab®
Regular Swab**
ITM-1148589 (white)
Minitip Swab
ITM-1066011 (blue)

LAB233 Anaerobic Culture
LAB240 Fungal Culture
LAB1231719 Fungal Stain
LAB235 Gonorrhea culture
LAB3417 Group B Strep, PCR
LAB3527 Group B Strep, Pen Allergy

LAB1231050 Group A Streptococcus PCR
LAB1231076 MRSA Screen PCR
LAB1231036 Respiratory Culture
LAB3691 Trichomonas Antigen (female)
LAB2111191 Wound Culture w/gram stain



**UTM®
Regular Swab**
ITM-1186125
Minitip Swab
ITM-1140078

LAB9520 Bordetella pertussis, PCR
LAB1230607 COVID-19 PCR (Batch)
LAB1231077 SARS-CoV-2, Flu A/B PCR
LAB1230746 COVID-19, Flu A&B, RSV PCR

LAB1231142 Flu A/B & RSV PCR
LAB9170 HSV PCR for Lesions
LAB9440 Mycoplasma/Ureaplasma Culture
LAB1230844 VZV PCR



OSOM® BVBlue®
Order via Lab
Client Services

LAB1231044 Vaginosis Screen



**24-Hour Urine
Jug****
ITM-1066912



Sterile Containers**
ITM-1072848

- ❖ Body Fluid ❖ Sputum
- ❖ Stool (unpreserved)
- ❖ Fresh Tissue ❖ Urine



Stool Kit **
ITM-1114530

Cary-Blair & SAF | Fill both
containers. | Always check
container label! Do not go
by color!



APTIMA®**

Unisex Swab Kit ITM-1065957 (white)
Multitest Kit ITM-1150249 (orange)
Urine ITM-1068485 (yellow)

- ✓ Patients 13 years and younger must be collected using the APTIMA® collection device. This testing will be sent out to a reference laboratory.
- ✓ Order Reference Miscellaneous testing and indicate test name in the order comment.
- ✓ Include source and specimen type.



Urine Tubes**

Tiger top UA Preservative Tube ITM-1118837

Urine C&S Kit Gray ITM-1065947

Yellow UA Urine Tube ITM-1065946




Occult Blood CRC Testing Vial**

polymedco FIT Kit

OC-Light® S FIT - Order via Client Services

Label must be filled out completely with patient name, date of birth and collection date of and time before mailing.

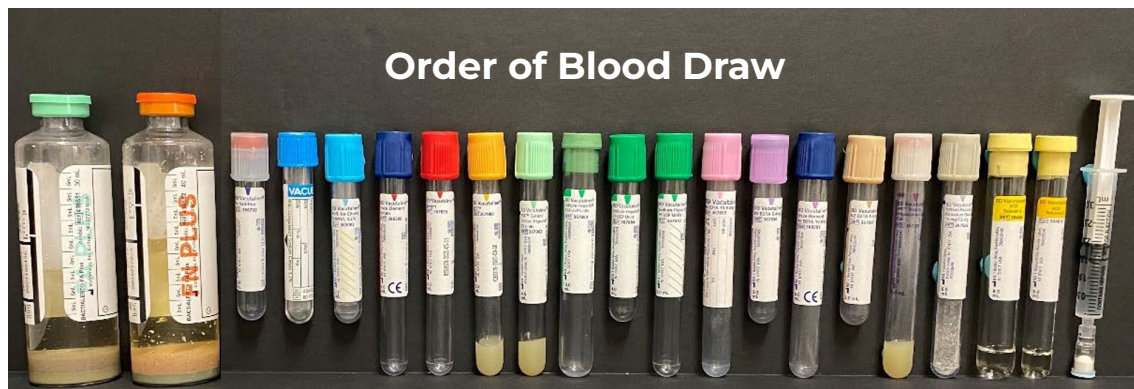

Order of Blood Draw

- ✓ Always review the lab test catalog** for important test collection information, including minimum volume, number of times the tube needs to be inverted after collection, if tube needs to be centrifuged/spun or left whole blood, specimen integrity times, and another collection and handling instructions.

Labeling

- ✓ **Every specimen must be labeled with patient's legal full name and date of birth.**
- ✓ Label should also include specimen type, date and time of collection, and source if needed.
- ✓ Do not place labels on the lid or bottom of the container, or stick labels to each other.
- ✓ Do not cover up patient info or barcodes on the labels.
- ✓ Do not cover up "windows" on tubes, leave colored edge of vacutainer visible.
- ✓ See lab catalog for more information on labeling.**

- ✓ Always use **2 Patient Identifiers** (full legal name and date of birth) on the specimen label and on the order.
- ✓ Always make sure that lids to containers, tubes and vials are on correctly and tight! Leaky specimens will be rejected.

GENERAL TEST LISTING, NOT ALL TESTS ARE LISTED**

****Review Lab Test Catalog: corewellhealth.testcatalog.org**